

KK-LAB



TOTAL WEIGHING SOLUTION™



Panasonic Kenko Electric. Co.Ltd



...Clearly a Better Value

Unit weights of up to 5 different samples can be easily entered, stored and recalled for use.



Comparator function

Compare samples to target values or pass/fail criteria and clearly indicate the results.



Formulation mode

Convenient for making many measurements of minute samples and seeking the total mass.



WindowsDirect Communication Function

Send balance data to Excel or other Windows applications without any data communication software installation required. By combining standard AutoPrint functions with typical spreadsheet functions, even difficult applications can be easily automated.

*RS232C interface is needed.



Panasonic Kenko Electric. Co.Ltd



JQA-0376



Capacity	82g	120g	220g	62g	120g	220g
Minimum Display	0.1mg					
Repeatability (Standard Deviation)	≤ 0.1mg					
Linearity	± 0.2mg					
Response Time*1	Approx. 3.0 seconds					
Operating Temperature and Humidity Limits	5-40°C 20~85%*2					
Temperature Coefficient for Sensitivity (10-30°C)	± 2 ppm/°C					
Pan Size (mm) approx.	φ 91					
Main Body Dimensions (mm) approx.	213(W) × 356(D) × 338(H)					
Main Body Weight (kg) approx.	6.2			6.0		
Power Requirement	12V, 1A					
Internal Calibration	●			—		

Capacity	Minimum Display		
300 gr	0.1 mg		
400 gr	0.1 mg		
1kg / 2 kg / 3 kg / 4 kg / 5 kg / 6 kg	0.01 & 0.001 gr		

*1 The response time is a representative v
*2 No condensation

Easy Setting Best fit to weighing application

Quickly adjust the desired ratio of stability and response for every application, even during measurement, with one-touch operation.



Option

Description
EP-80 Printer
EP-90 Printer
RS-232C interface IFB-102A-UNC
RS-232C Cable 25P-9P
Protection Cover
STABLO-EX
USB conversion kit

Expanded Piece Counting function

Unit weights of up to 5 different samples can be easily entered, stored and recalled for use.

Comparator function

Compare samples to target values or pass/fail criteria and clearly indicate the results.

Formulation mode

Convenient for making many measurements of minute samples and seeking the total mass.


WindowsDirect Communication Function

Send balance data to Excel or other Windows applications without any data communication software installation required. By combining standard AutoPrint functions with typical spreadsheet functions, even difficult applications can be easily automated.

*RS232C interface is needed.



Dimensions

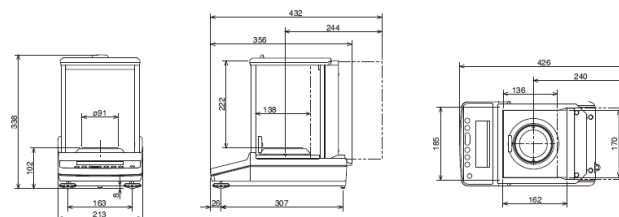

Panasonic
KENKO®



...Clearly a Better Value



LIT: FG-K 01/08

Panasonic Kenko Electric Co, Ltd
36 Sisenju, Okio, Osaka 234-23354 Japan

Tel. 81-32-565434667

Specifications subject to change without notice

